3H *Smart*

Contents

3H LED MASK Premium

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
Contents

3H LED MASK Premium

1. 3H LED Mask Accessories & Effects	12
2. User Instruction/How to use	13
2-1. How to use LED Mask	
2-2. How to use the Controller	
※ LED Mask Clinical Trial Result	
3. Product Precautions	16
- Product Warnings and Cautions	
- In-Use Injury Warning	
- Warnings about using the product	
- Warnings and Cautions for Adapters	
4. Troubleshooting and Q&A	18
5. Product Specification	19

3H LED MASK Premium

1. 3H LED Mask Accessories & Effects

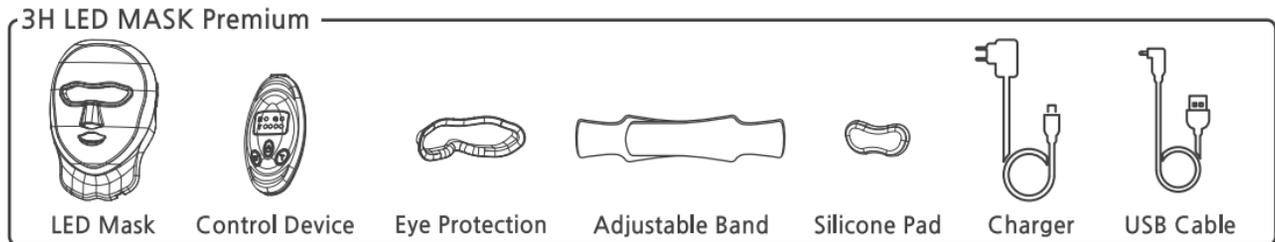Red, Near-Infrared Ray, Three-wavelength LED of Derma LED

- 3H LED Mask has LEDs with 3 wavelengths, which helps to improve facial skin elasticity, skin tone and tighten dermal density.
- The complex wavelength of LED (Red LED)/Near-Infrared Ray LED (N-IR LED) / Combi LED (Combi LED) are transmitted from the skin surface to the dermal layer

Eye protection design by 3H LED Mask.

- Passed & Certified with the International Photobiological Safety Test Standard when directly exposed to LED light.
- EC/EN 62471 International Standard Certification by Intertek.

Wear comfortably & Use conveniently.

- Passed & Certified with the International Photobiological Safety Test Standard when directly exposed to LED light.
- EC/EN 62471 International Standard Certification by Intertek.

2. User instruction/How to use

2-1. How to use LED Mask

1. Fully charge the controller before use. The charge status of the battery is indicated by a 4 step indicator.
2. Before use, wash the face thoroughly with facial cleanser and apply the included Multi-peptide Cream.
3. Connect the strap band to the LED Mask.
4. Connect the cable from the power terminal of the LED Mask to the USB terminal on the controller.
5. Adjust the LED Mask to the to opening positions of eyes, nose, mouth, then adjust the strap band.
6. Press the power button on the controller for about 2 seconds to turn on the LED Dmask.

Then LED Mask with light up with a beeping sound.

7. Press the time setting button briefly to set the operation time.


(Can be set up to 20 minutes and 5 minutes of interval.)

8. After setting the timer, press the power/operation button briefly to activate. (The LED will light up.)
9. At the end of the preconfigured operation time, the LED Mask will automatically will power off.

* Recommended use time per day: 10 to 20 minutes per day.

If exceeding daily use of the recommended time may cause the skin to feel warm & may cause redness.

2-2. How to use the controller

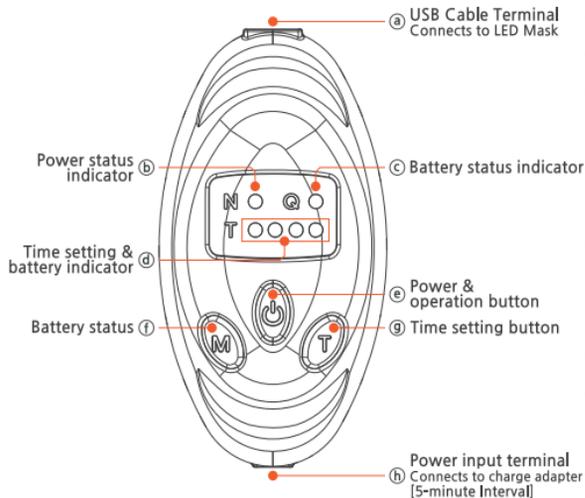a USB cable terminal : Led Mask connection cable terminal

b Power status indicator : Lights in blue when the controller is operating

c Charging status indicator : Lights in red when the charging adapter is connected to the power input terminal.[The indicator turns off automatically when the battery is fully charged.]

d Battery status : Press the M button once and it blinks three times and the remaining battery status is displayed.

e Power&operation button : Press&hold the button for about 2 seconds to power on. After setting the time/mode, a short press of the button will activate the LED Mask. [the operation of LED Mask will end and the controller automatically shuts off]

f Time setting button : Each press of the T button, it sets to 5~20minutes in 5minute increments.

g Power Input Terminal : Connect the charging adapters to charge the controller battery.

- The controller is automatically charged when connected to the charger, The LED Mask operates during charging and is automatically powered off to prevent battery discharge after charging is complete.

- When fully charged, all 4 indicators be shown and lit.

- LED Mask charges and operates when connected to 5V/1A portable battery.

3H LED MASK Premium Clinical Trial Result

P&K Skin Research Center of Korea. Global Frontiers of Skin Research.

0% Side Effects
Proven Effectiveness of
Skin Improvement


Skin
Moisturizing
96.5%

Smoother Skin
89.3%

Skin Texture
Improvement
98%


Skin Elasticity
96%

Skin Elasticity

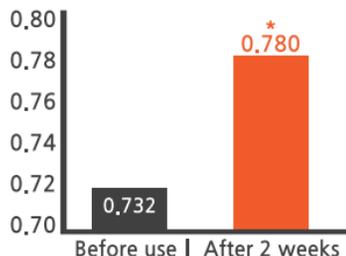
Skin elasticity increased by 6.631%
Purpose: Improving Skin Elasticity

Skin Moisturizing

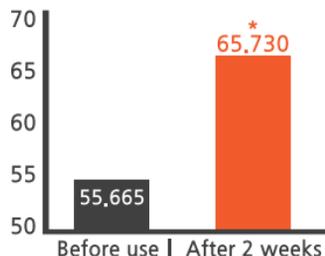Skin moisturization increased by 18.082%
Purpose: Improving Skin Elasticity

Dermal Density

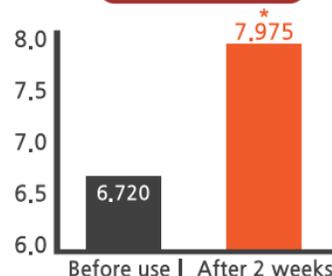Dermal Density increased by 18.397%
Purpose: Improving Skin Elasticity

Use: 20 mins daily use for 2 weeks / Subject: 22 of 30 to 50 year old women /
Evaluation: Test for skin elasticity, skin moisturization, improvement of dermal density

3H DIY Home Near-Infrared LED Mask Premium

Anti-aging

Skin Tone

Skin Elasticity

Anti-wrinkle

Anti-blemish

Advanced Esthetic Management. Now Experience at home.

3. Product Precautions

Product warnings and cautions

- Do not aim the LED light and avoid looking directly at the LED light as it may cause temporary discomfort.
- Intended for adults over the age of 18. Not recommended for infants or children.
- If you are pregnant or suspect you are pregnant, please consult your doctor before use.
- This product contains a polymer battery. Please consult your local requirements regarding disposal of the controller.
- If the product heats up or becomes uncomfortably hot, discontinue use.
- When charging is complete, please disconnect from the outlet.
- Do not store in direct sunlight or at temperatures above 100 degrees Fahrenheit (38 degrees Celcius).
- Use only the supplied universal USB cable and the adapter. Using incorrect power may cause damage to the product.
- Do not open or repair the LED Mask or controller.
- Dropping the product on a hard surface may cause damage to the product. Damaged products may cause harm to the skin.

Product usage warnings and cautions

- Handle the power plug with a dry hand. May cause electric shock.
- Be sure to hold the end of the cable, not the cable wire, when pulling out the cable.
- Can cause static electricity when used in dry environments.
- If contaminated, do not wash the product with water, but wipe it with a cloth. If the contamination is serious, use a little detergent.
- If the product is wiped with alkaline cleaner, thinner or alcohol, it may damage the product.
- If not used for an extended period, charge at least once every three months to maintain battery life.

Warnings about using the product

- If you are in the process of dermatologic treatment, or if your skin condition is abnormal, please consult with a dermatologist before using the product.
- If you have allergies or have a sensitive skin, please use after consultation with a dermatologist.
- Do not use this product if you have a skin condition or irritation, such as a sore skin or severe acne condition, burns, or bacterial skin infection.
- When using the product, depending on the individual skin condition, the skin may be irritated, such as reddening or tingling.
- If such irritation persists, discontinue use immediately and consult a dermatologist.
- Do not use this device if you are in the skin recovery period due to skin-related surgery, plastic surgery, or hair loss procedure.
- This product is a face-only product and is not recommended for use on other parts of the body.
- During Pregnancy, menstruation, lactating women are not recommended for use. Temporary skin troubles may occur.
- Do not use the product for any purpose other than what is described in this user's guide.

Warnings and cautions for adapters

- Use only the indicated power and voltage. If the connection is unstable, it may cause fire.
- If you notice a smoke or smell unpleasant odor from the product, discontinue use and contact Customer Service. Continued use may cause fire or electric shock.
- Do not damage, process, bend, twist, or heat the charging adapter and USB cable. Damage to the cord may cause fire or electric shock.
- Do not insert metal object or other conductors on the other side of the cord while the charging adapter is connected to the outlet.
- When connecting the controller and adapter, check the cable terminal direction. It may cause malfunction of the product.
- When unplugging the charging adapter from the wall outlet, grasp the adapter body and pull out. If pulled by the cable, it may cause damage or disconnection.
- Do not immerse in water or wash the charging device. It may cause fire or electric shock.
- Do not use a third-party charger other than the supplied charging adapter.

4. Troubleshooting and Q&A

Warranty

The warranty period for this product is one year after purchase.(However, the defect caused by consumer negligence is excluded.)

Warranty

The warranty period is calculated from the date of purchase, and if not confirmed, the warranty term will be applied based on the two months from the date of manufacture.

The following cases are not covered even within the warranty period.

- Foreign substances in the product (liquids such as water and beverages)
- For products that are not electronically suitable for specified countries
- For products distributed through unauthorized channel
- Damage and failure due to natural disaster (lightning, fire, flood, earthquake)
- Physical damage and malfunction due to consumer negligence
- Damaged using other non-compatible and/or third-party products other than the dedicated components

Warranty claim checklist

- If the product doesn't power on : The battery may be drained. Use after charging.
- If malfunctions : Check for foreign objects/materials on the connection terminal and reconnect.
Make sure the power button is pressed for two seconds.
- When the battery does not charge : Replace the charging adapter to charge.
- If the battery life is different from the user manual : Lifetime may be shortened when not operated for a long time.
Because it is a consumable item,
it may need to be replaced when performance degrades.

Model Number		1-year Limited Warranty for product
Purchase Date		6-months Limited Warranty for battery

5. Product specification

Power specification	Input AC 110 ~ 240V 60Hz, Output DC 5V 1.2A
Battery	3.7V 1,000 mAh
Charging time	Approximately 1 hour
Number of use per full charge	Approximately 4 times per charge (based on 20 minutes of usage)
LED wavelength	3H LED MASK Premium 3 wavelengths
Time settings	5 minutes, 10 minutes, 15 minutes, 20 minutes
Dimension	7.87" x 11" x 3.54" / 200mm x 280mm x 90mm
Weight	LED Mask: 280g
Manufacturer	Teddy Korea. Co., Ltd. 312-1, ACE hitech city 2, Mullaedong 6 ga 25, Seonyu-ro 13-gil, Yeongdeungpo, Seoul, Korea
Distributor	3H Co.,Ltd 140-10(Yuram-dong) Yuram-ro, Dong-gu, Daegu Metropolitan city, Korea
Customer Service	+82-80-575-3333

Direct Power Supply
HU10772-17008
MSIP-REM-YH6-SPE-TU120
Battery Safety Certification
YU10770-18001B
Electronic Wave Certification
R-R-TDK-TKB-FG03

Photobiological Safety Test
IEC/EN 62471
170900067SEL-001

*Product specifications and specification designs are subject to change without prior notice.

*This product can be compensated for by the insurance of the goods. (K-Biz Small Business mid-time)

3H Co.,Ltd.

www.3hk.co.kr

MADE IN KOREA